

CAPOVARIO HELMETS

VERSATILE HEAD
PROTECTION




SUPPORTS & BRACES HELMETS

Just like a fingerprint, every head is unique—so the Capovario helmets are designed with maximum customisation in mind. The adjustable head circumference, crown padding, and Y-straps ensure a precise and comfortable fit.

Customise your Capovario helmets with a range of accessories, a chin guard, a closed top, or buckle and strap options. Available in five sizes and multiple colours, it features a durable, anatomically shaped plastic base that adapts to your head for a secure fit.

Crafted from high-quality microfiber fabric, this helmet is breathable, washable at 40°C, and can be disinfected. The removable inner padding can also be easily washed or replaced as needed.

- Continuously adjustable head circumference for a perfect fit
- Removable inner padding for easy cleaning and customisation
- Washable at 40°C and can be disinfected for optimal hygiene
- Compatible with innovative accessories like chin guard, visor, and more

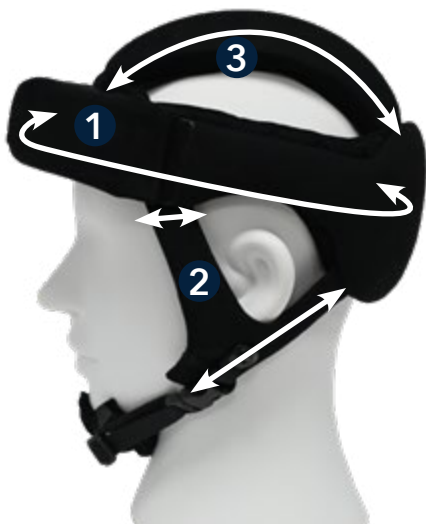

Areas that can be adjusted

1. Head Circumference
2. Position and depth of Y-Straps
3. Depth of crown pads



'Customisable Head Protection for a Perfect Fit'



SIZING CAPOVARIO

Getting the right size for your Capovario helmet is essential for both safety and comfort. A properly fitted helmet ensures optimal protection by staying securely in place during impact, while also preventing pressure points or discomfort during wear.

Capovario helmets offer a range of sizing options, and their customisation features—such as adjustable padding systems and precision fit dials—allow wearers to fine-tune the fit for their individual head shape.

These customisable elements make it easier to achieve a secure, personalised fit, even if you're between standard sizes. Taking the time to get your sizing right, with the help of these adjustments, enhances both performance and protection.



CAPOVARIO HELMETS

The Capovario helmet offers adaptable, comfortable head protection for patients with a wide range of medical needs. Standard integrated forehead protector helps shield the face and nose.

Available in five sizes, the Capovario ensures optimal protection across different age groups and head shapes. Patients can also choose from a range of attractive colours.

Crafted for hygiene and comfort, the Capovario features a high-quality microfibre cover that is breathable, washable at 40°C, and easy to disinfect. Each standard helmet includes the helmet shell, Y-straps, and one set of inner padding.

CAPOVARIO ACCESSORIES

To further enhance the functionality and patient comfort of the Capovario helmet, MedFac offers a range of practical accessories designed to meet individual care needs.

Buckle Cover

This soft, protective cover helps cushion the fastening buckle, reducing pressure and skin irritation—ideal for patients with sensitive skin or prolonged helmet use.

Chin Strap

The adjustable chin strap ensures a secure fit, providing additional stability and peace of mind, especially during active daily routines or therapeutic activities.


Closed Top

For patients who require complete cranial coverage, the closed top accessory adds a protective layer over the head, offering extra safety without compromising comfort.

Inner Padding

Replacement inner padding is available to maintain hygiene and extend the life of the helmet. The padding is breathable, removable, and washable, ensuring comfort and cleanliness over time.

All accessories are designed with the same focus on safety, hygiene, and adaptability as the helmet itself, allowing care teams to tailor the Capovario to the specific needs of each patient.

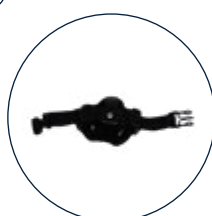

CHIN STRAP



INNER PADDING



CLOSED TOP



BUCKLE COVER

CAPOVARIO HELMETS & ACCESSORIES

PRODUCT CODE	DESCRIPTION
H600500-P	Capovario Size 0. Pink
H600500-T	Capovario Size 0. Turquoise
H600501-P	Capovario Size 1. Pink
H600501-T	Capovario Size 1. Turquoise
H600502-S	Capovario Size 2. Black
H600503-S	Capovario Size 3. Black
H600504-S	Capovario Size 4. Black
H600515-1	Chin protector Size 0,1 & 2
H600515-2	Chin protector Size 3 & 4
H600525-2	Closed top Size 2
H600525-3	Closed top Size 3
H600525-4	Closed top Size 4
H600540-0	Inner Padding Size 0
H600540-1	Inner Padding Size 1
H600540-2	Inner Padding Size 2
H600540-3	Inner Padding Size 3
H600540-4	Inner Padding Size 4
H600580-2	Buckle Cover

CAPOVARIO-SOFT HELMETS

The Capovario-Soft helmet is designed to protect patients from head injuries caused by uncontrolled or unsafe movements. Made from durable cotton fabric and padded with special integral foam, it also includes an impact-resistant foil layer for added safety.

Available in five adjustable sizes, the helmet allows for a head circumference adjustment of approximately 4-5 cm and can be further customised for vertex height. Velcro®-secured forehead and neck protectors offer added flexibility, while the improved Y-strap system ensures a secure, comfortable fit tailored to each individual.



CAPOVARIO-SOFT ACCESSORIES

To enhance comfort and safety, the Capovario-Soft helmet can be paired with a range of practical accessories tailored to individual needs:

Buckle Cover

Softly padded to reduce pressure and skin irritation around the fastening area, offering extra comfort during prolonged wear.

Chin Protector

Provides gentle cushioning for the chin strap area, helping to prevent rubbing and ensuring a secure, comfortable fit.

Neck Protector

Easily attached via Velcro®, this accessory adds extra coverage and support to the back of the head and upper neck, increasing protection during backward movements or falls.

Forehead Protector

Designed for patients who need additional frontal protection, this cushioned guard helps shield the forehead while maintaining breathability and flexibility.



BUCKLE COVER



CAPOVARIO-SOFT HELMETS & ACCESSORIES	
PRODUCT CODE	DESCRIPTION
H600100-02	Capovario-Soft Size 0 Navy
H600100-12	Capovario-Soft Size 1 Navy
H600100-22	Capovario-Soft Size 2 Navy
H600100-32	Capovario-Soft Size 3 Navy
H600100-42	Capovario-Soft Size 4 Navy
H600000-62	Capovario-Soft Chin Protector Navy
H600000-72	Capovario-Soft Forehead Protector Navy
H600000-82	Capovario-Soft Neck Protector Navy
H600580-1	Buckle Cover

MEASURING FOR HELMETS

The Capovario helmet features an anatomically shaped shell made from sturdy plastic, allowing the head circumference to be adjusted by up to 4 cm. This flexibility helps achieve a secure and personalised fit for each patient.

To select the right size:

- Measure the head circumference just above the eyebrows and around the widest part of the head using a soft tape measure.
- Refer to the size charts below to choose the size range that includes the measured circumference.
- Allow for adjustability, as each helmet size offers up to 4 cm of fine-tuning for a snug and comfortable fit.
- A properly fitted helmet ensures maximum protection and optimal wearability throughout the day.



CAPOVARIO			
SIZE	HEAD CIRCUMFERENCE	WEIGHT	COLOURS
0	44 - 48 cm	270 g	Pink/ Turquoise
1	48 - 51 cm	280 g	
2	52 - 56 cm	290 g	Black
3	56 - 60 cm	305 g	
4	58 - 62 cm	320 g	

CAPOVARIO SOFT			
SIZE	HEAD CIRCUMFERENCE	WEIGHT	COLOURS
0	43 - 47 cm	110 g	Navy
1	47 - 49 cm	120 g	
2	49 - 52 cm	130 g	
3	52 - 56 cm	145 g	
4	56 - 60 cm	150 g	

ADJUSTING YOUR HELMET

When delivered, the helmet is set to its largest size. Follow these steps to adjust it for a personalised fit:

Position the Helmet

Place the helmet on the head. It should sit level, about a finger's width above the eyebrows, with the ears exposed.

Adjust Head Circumference

Loosen the two thumbscrews under the Velcro forehead padding. Slide the plastic tabs evenly on both sides for a snug fit. Tighten the screws to secure.

Check the Fit

Rotate the helmet gently with both hands. If the skin on the forehead moves with it, the fit is correct. If the helmet shifts independently, repeat step 2.


Adjust the Straps

Open the Y-strap adjusters and position the straps securely. Fasten the chin buckle, leaving no more than one finger's width of space. Adjust the padded strap with the Velcro underneath as needed.

Final Fit Check

Pull the helmet forward and backward. If it blocks vision or slips, re-tighten the straps.

For further assistance, contact your medical supply provider.

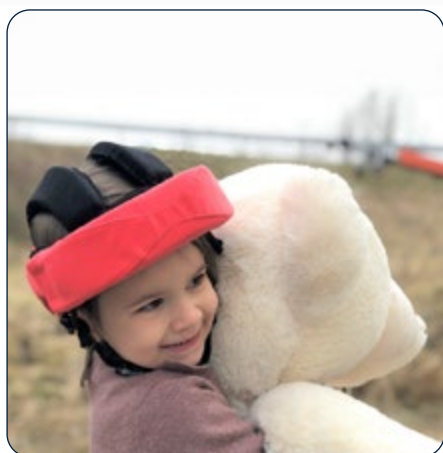

Visit us online
www.kinetecuk.com

Or contact us directly for your Orthotic requirements.

UK Enquiries:

MedFac UK Ltd
7 Redan Hill Industrial Estate
Redan Road
Aldershot
GU12 4SJ

Tel: +44 3330 347400
Fax: 01252 329415

enquiries@kinetecuk.com
www.kinetecuk.com

International Enquiries:

Kinetec SAS
Zone Industrielle de Tournes
Rue Maurice Périn
F-08090 Tournes
France

Tel: +33 324 29 85 05
Fax: +33 324 33 51 05

contact@kinetec.fr
www.kinetec.fr

MedFac
ORTHOTICS