

Unique
veterinary
casting & splinting
materials

 **Vet-LITE™**

by  runlite



Vet-Lite™ Classic Splinting

Rolls & Sheets



Vet-Lite Classic is a strong material that can be used for any cast or splint of your choice. It stretches and conforms well to contours for a perfect fit.
Comes in sheets and rolls.

VET-LITE CLASSIC ROLLS	
902	Vet-Lite Roll 5 x 180 cm (20/Pkg)
903	Vet-Lite Roll 7,5 x 180 cm (20/Pkg)
904	Vet-Lite Roll 10 x 180 cm (10/Pkg)
905	Vet-Lite Roll 15 x 180 cm (10/Pkg)
VET-LITE CLASSIC SHEETS	
993	Vet-Lite Sheet 7.5 x 37,5 cm (15/Pkg)
994	Vet-Lite Sheet 10 x 37,5 cm (15/Pkg)
996	Vet-Lite Sheet 15 x 37,5 cm (15/Pkg)

Vet-Lite 3D™ Splinting

Sheets



Vet-Lite 3D material features a base with 3-Dimensional Polyester fabric which contributes to increased strength. Available in sheets in three different thicknesses - 3mm, 5mm and 7mm.

VET-LITE 3D SHEETS	
9896	Vet-Lite-3D Sheet 600 x 600 x 3mm, Beige
9897	Vet-Lite-3D Sheet 600 x 900 x 3mm, Beige
9898	Vet-Lite-3D Sheet 600 x 600 x 5mm, Beige
9899	Vet-Lite-3D Sheet 600 x 900 x 5mm, Beige
9800	Vet-Lite-3D Sheet 600 x 600 x 7mm, Beige
9801	Vet-Lite-3D Sheet 600 x 900 x 7mm, Beige

Manosplint® - Montana & Ohio-T

Thermoplastic sheets

Montana

MONTANA				
THICKNESS	PERFORATION	COLOUR	SHEET SIZE	
			DEMI 12" x 18" (30.5 x 46cm)	STANDARD 18" x 24" (46 x 61cm)
1/8" (3.2mm)	Solid	White	M1MN4D00WH	M1MN4S00WH
3/32" (2.4mm)	Solid		M1MN3D00WH	M1MN3S00WH



Memory

Moderate. Will spring back slightly when heated.

Resistance to Stretch (RTS)
 Moderate to maximum resistance to stretch. Can be stretched and pulled around curves without thinning or losing control.

Drape

Moderate. Takes on the contours of the target area for an intimate fit by using only moderate pressure.

Rigidity
 Maximum. Splints will hold their shape against the pull of strong dynamometers and hypertonicity.

Bonding & Finish

Uncoated. Slightly sticky surface when warm. A permanent bond can be easily formed by applying dry heat and pinching it together.

Surface Finish & Edging
 Resists fingerprinting and edges beautifully. Leaves sealed, smooth edges when cut warm.

Ohio-T



OHIO T			
THICKNESS	PERFORATION	COLOUR	STANDARD 18" x 24" (46 x 61cm)
1/8" (3.2mm)	Solid	Black	M1OT4S00BK
3/32" (2.4mm)	Solid	Black	M1OT3S00BK
1/16" (1.6mm)	Solid	Black	M1OT1S00BK

Memory

100% memory - when fully reheated can regain original size and shape.

Resistance to Stretch (RTS)
 Moderate handling with smooth strokes during moulding. Shapes easily.

Drape

Moderate pressure needed for a conforming fit.

Rigidity
 Maximum - strong support.

Bonding & Finish

Semi-coated material allows for it to 'pinch and pop'. Acts as a second pair of hands to keep product in place during application. Resists fingerprints with a shiny finish.

Hook, Loop and Strapping Solutions




Vet-Lite Velcro T/Hook
An easy solution to attach your pile closure straps to your splint or brace. Bonded to the back of the 'hook' is the same cotton gauze used in the Finishing/ Edging Tape which bonds instantly and permanently on all our materials
Product code
641970000 Black, 5 x 50 cm (5 sheets)



Low profile hook
Self-adhesive hook rectangle pieces. Works with all loop straps. Smooth to the touch. 80 x 38mm. 100/pkg.
Product code
M600000007

PRO RANGE	
M600000042	Manosplint Pro Adhesive Hook 2.5cm x 25m Black
M600000043	Manosplint Pro Adhesive Hook 2.5cm x 25m White
M600000044	Manosplint Pro Adhesive Hook 5cm x 25m Black
M600000045	Manosplint Pro Adhesive Hook 5cm x 25m White
M600000046	Manosplint Pro Sew On Loop 2.5cm x 25m Black
M600000047	Manosplint Pro Sew On Loop 2.5cm x 25m White
ESSENTIALS RANGE	
M600000048	Manosplint Essentials Sew On Loop 2.5cm x 25m Black
M600000049	Manosplint Essentials Sew On Loop 2.5cm x 25m White
M600000050	Manosplint Essentials Sew On Loop 5cm x 25m Black
M600000051	Manosplint Essentials Sew On Loop 5cm x 25m White

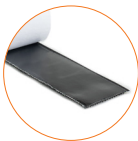

Mikrofix
Self-adhesive Hook and Loop. Ultra-thin, strong and high quality. The Hook material has very short hooks which does not tear unnecessarily and gives the patient greater comfort. Supplied in 50mm width and 25m length.
25041 0000 Mikrofix, hook (black)



Black Elasticated Loop
Stretchy loop. Knitted. Cuts easily. Smooth to the touch. Works with all hook straps. 1" x 41' (25mm x 12.5m)
Product code
M600000028 25mm x 12.5m



Manosplint® Soft Foam Loop
Compatible with hook tapes, non-stretch strapping material made with soft, comfortable covering over foam core. Thickness 3mm
Dimensions: 5cmx9m
Product code
M600000063 Black



Strap Elastic
Material provides both stretch and padding. Hook receptive on both sides. Works great with our double-sided hook dots!
Maximum comfort and moderate elasticity. Beige only. 10/pkg.
Product code
M600000055 38 x 295 mm

Heating Solutions

HEAT GUN

Outstanding, powerful 2000W heat gun.
Designed to be used both upright and hand-held.

Features include

- Wide, flat base provides superb stability when placed upright
- Easy-access, big button control even when in upright position
- Two heat and airflow settings for preferred heating methods
- Lightweight, ergonomically designed for optimal weight distribution
- Range of nozzles provided
- Two year warranty

Product code **M203515181**



PORTABLE ELECTRIC SPLINT PAN (30CM SQUARE)

The pan simply plugs into the mains to deliver 1500W of power at temperatures of up to 240°C. It's a versatile, controllable pan with an inbuilt thermostat that enables you to control output.

Features include

- Cool-touch handles and lid knob. Powerful 1500W
- Durable and strong aluminium body with glass lid and steam vent
- Non-stick and large surface area: 30cm x 30cm x 5cm
- Inner dimensions: 30cm x 30cm x 5cm
- External dimensions: 32cm x 45cm x 13cm

Product code **M202007089**




BI-FOLD LID HEAT PANS



5 Litre

Bi-Fold Lid **M200000054**

The Bi-fold lid heat pans are designed to heat thermoplastic splinting sheets in warm water while maintaining a controlled temperature. This allows to soften the material safely and efficiently so it can be quickly moulded into custom splints.

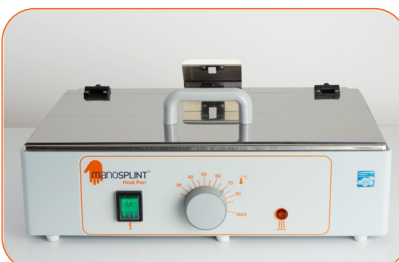

9 Litre

Bi-Fold Lid **M200000055**

Features include:

- Stainless steel construction
- Thermostatic temperature control
- Folding 'bi-fold' lid that saves space
- Insulated handles
- Perforated tray inside the bath

Cutting Solutions



Utility Knife

- Ergonomic grip for easier cutting and firm control
- Retractable blade
- Easy to change blades

Product code

M400000022

Replacement Blades

M400000023 10/pkg



Steel Scissors (Small)

- Pointed, blades for perfect cutting of most materials
- Ideal for cutting strapping materials as can be placed in the water bath for cleaning
- Simple steel handles
- Total length 15cm; blade length 6.5cm

Product code **M40112015**



Classic Universal Scissors

- Cuts a range of materials easily
- Large, ergonomic handles for easy gripping and control
- Strengthened plastic handles for robust use
- Total length 20cm; blade length 8.5cm

Product code **M40762220**



Super Cast Scissors

High quality cast scissors which easily cut through all materials. The lower blade has fine serrations and a sharp-angled upper blade. Stainless steel and autoclavable.

Product code 18cm

512070000

Product code 14cm

512080000



Steel Shears

- Robust and easy to clean
- Supplied either straight or curved
- Length 21cm; blade length 7.5cm

Product code Straight Blade

M400000013

Product code Curved Blade

M400000012



Heavy Duty Shears

- Ideal for cutting cold plastics up to 3.2mm thick
- Self-opening
- Short, sturdy blades hot forged for endurance
- Length 24cm; blade length 4.5cm

Product code **M40113824**

Stockinettes



Cotton/Poly Stockinette

White only

25 yard

(22.86m) rolls

M500000022 2" (5cm)

M500000023 3" (7.6cm)

M500000024 4" (10cm)



Cotton Stockinette

Black

25 yard (22.86m) rolls

M500000025 2" (5cm)

M500000026 3" (7.6cm)

M500000027 4" (10cm)

Padding & Edging



X-Lite® Foam

- Padding is self-adhesive.
- Supplied in white and 2.75m length.

64480 0003 (70mm, 20 rolls)



ViscoFoam Padding

Self-adhesive, slow recovery foam padding in

3/8" x 16" x 24" (9.5mm x 41cm x 61cm).

M300000001 Medium Density, Blue

M300000002 Soft Density, Pink



CushionFoam Padding

Self-adhesive, open-cell, soft padding with 4-way stretch. 46cm x 61cm. 4 sheets pack. Choose between 3.2mm or 6.4mm thickness.

M300000027 3.2mm, Black

M300000030 6.4mm, Black



X-Lite® Adhesive Foam Edging

A soft and comfortable adhesive foam edging for X-Lite splints.


64238 0000 25mm x 4m (Black)



Vet-Lite Plus Edging

Vet-Lite Plus Edging tape gives a neat finish to the splint and secures a soft edge. Supplied in a dispenser box.

64238 0000 25mm x 4m (Black)

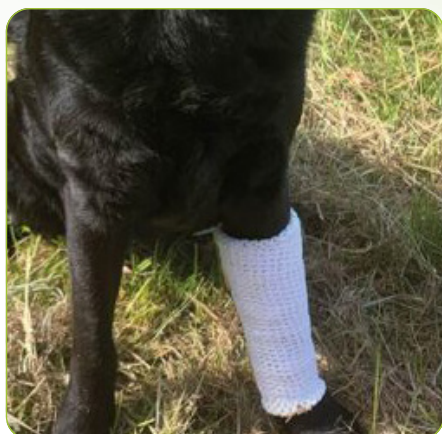

Visit us online
www.runlite.com

UK Enquiries:

MedFac UK Ltd
7 Redan Hill Industrial Estate
Redan Road
Aldershot
GU12 4SJ
Tel: +44 (0) 3330 347400

vetlite@runlite.com
www.runlite.com

 **Vet-LITE™**

brought to you by

 **runlite^{SA}**